



High sensibility and accuracy with TRFIA



Test kits transport at 4-30°C



Small sample requirements



MNCHIP - Immunology

Immunoassay Analyzer Pointcare® iM1

Features & Benefits

Easy to use

Only 5 operation steps
Compact design and light weight

Express detection

Deliver results in 5-10 minutes

Reliable accuracy

Time-resolved fluoroimmunoassay (TRFIA)
Integrated quality control system (RQC)


Cost-effective

Test kits can be stored at 4-30°C, saving shipping cost
Affordable analyzers and test kits


Operation




Add 2µl sample into the diluent tube

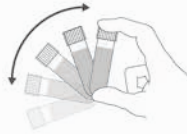

Gently invert the collection tube several times just prior to sample transfer

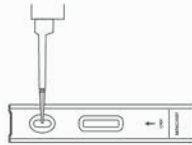

Transfer 100µl diluted sample to the sample hole

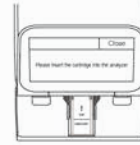

Insert the cartridge into the device



Read result

*Remarks: the operation may differ according to different types of samples or reagents. Please refer to package inserts for specific operating procedures.

Test Kits

Inflammation	PCT Test Kit
	CRP Test Kit
	SAA Test Kit
Cardiovascular	cTnl Test Kit
	Myo Test Kit
	BNP Test Kit
	H-FABP Test Kit
Coagulation	D-Dimer Test Kit
Diabetes & Renal	HbA1c Test Kit
	CysC Test Kit
	β2-MG Test Kit
	NGAL Test Kit
...	

* More test kits coming soon...

Specifications

Sample Types	Serum or plasma * Please refer to package inserts for details.
Interfaces	USB
Memory Capacity	Up to 50,000 test results
Power Bank	Capable
Dimensions	Dimensions: 210mm(L)×125mm(W)×175mm(H) Weight: 2.0kg
Environmental Specifications	Operation Temp: 10-30°C
	Operation Humidity: 40-85%
	Power Supply: 24V