



Immunofluorescence



Test kits transport
at 4-30°C



Small sample requirements



MNCHIP - Immunology

Immunoassay Analyzer Pointcare® Reader V2

Features & Benefits

Easy to use

- Small in size and applicable to multiple scenarios
- a 4.3-inch screen for intelligent and convenient operation

Express detection

- Deliver results in 5-15 minutes.
- Under the "fast detection" mode, results are reported immediately, allowing for consecutive testing

Reliable accuracy

- Advanced immunofluorescence quantification technology provides high-quality quantitative detection
- Integrated quality control system (RQC)


Cost-effective

- Test kits can be stored at 4-30°C, saving shipping cost Affordable analyzers and test kits


Operation




Add sample into the diluent tube

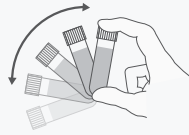


Gently invert the collection tube several times

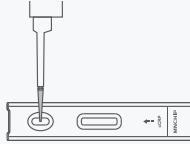

Transfer 100µl diluted sample to the sample hole

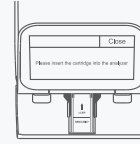

Insert the cartridge into the device

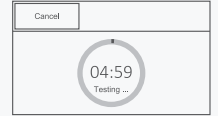

Read result

*Remarks: the operation may differ according to different types of samples or reagents. Please refer to package inserts for specific operating procedures.

Test Kits

Inflammation	cCRP
	fSAA
Fertility	cProg
Renal	cCys C
	fCys C
	cSDMA
	fSDMA
Cardiovascular	fNT-proBNP
	cNT-proBNP
	cTnI
Pancreas	cPL
	fPL
Endocrine	cCortisol
	cTT4
	cTSH
	HbA1c
Coagulation	cD-Dimer
...	

More test kits coming soon...

Specifications

Sample Types	Serum, plasma or whole blood. * Please refer to package inserts for details.
Time to Results	About 5-15 minutes
Connection Interfaces	WLAN, USB, Ethernet interface, compatible with LIS
Print Mode	Optional external thermal printer, office printer (optional), MDMP, LIS
Memory Capacity	Up to 50,000 test results
User Interface	4.3 Inch touch screen
Dimensions	Dimensions: 210(L)×125(W)×175(H)mm
	Weight: About 2.0kg
Environmental Specifications	Operation Temp: 10-30°C
	Operation Humidity: 40-85%
	Power Supply: 24V

Corporate Headquarters

MNCHIP Technologies Co.,Ltd
1-4F,Area,No.122 Dongting Rd,
Development Zone,300457 Tianjin P.R.China

Contact Us

sales@mnchip.com
service@mnchip.com
www.mnchip.com

Ver 2.1

MNCHIP All rights reserved.