



Xenon light source



Microfluidic Reagent Rotor

12 minutes

Sample-to-answer in 3 simple steps
and approximately 12 minutes

Clinical Chemistry⁺ – Celercare™ M5 On-site Blood Chemistry Analyzer


EASY TO USE

Fully automatic system – no special operating skills required


- 8.0-inch capacitive touchscreen with a user-friendly interface.
- Requires only 100 µL of lithium-heparinized whole blood, plasma, or serum.
- Barcode-based calibration information.
- Built-in printer for direct result printing.
- LIS compatible, no need to manually enter patient information.
- Results can be printed through an office printer or through data management platform (MDMP).

QUICK RESULTS

Complete test results in 3 simple steps, in approximately 12 minutes.

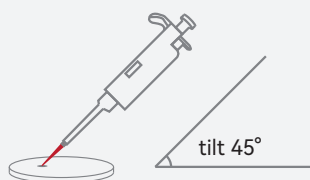

1. Add sample

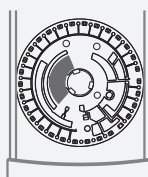

2. Insert disc

	Page Up	Page Down	Open	
Name: Age:				
ID:		Gender:		
Department:		Patient Area:		
Sample type:				
Lot:				
Item Name	Result	Indicator	Ranges	Unit
ALB	45.0		40-50	g/L

3. Read results

ACCURATE AND RELIABLE

Advanced technology ensures precise results


- Microfluidic discs with pre-installed reagent pellets ensure accurate analysis of blood samples and reagents.
- Stable measurement optics include a stroboscopic xenon lamp, a wavelength selection system, and a multiple-wavelength detector.
- Integrated quality control software monitors the entire process in real time (ensuring consistent performance of the analyzer, discs, and reagents, as well as accurate analysis of blood samples).

Specifications

Sample Types	lithium-heparinized whole blood, plasma, or serum
Sample Volume	100 µL
Detection time	Approximately 12 minutes
User Interface	8.0-inch touch screen
Reaction Temperature	37°C±0.2°C
Power Supply	120VA
Dimensions	260(L)×230(W)×325(H)mm
Net weight	Approx. 5.5kg

Parameters	Up to 19 parameters in one test, 33 parameters configured into 10 profiles
Reagent Disc	Disposable, prepackaged with self-contained lyophilized reagent
Calibration	Automatic self-calibration by scanning QR Code on the aluminium foil pouch
Quality Control	Built-in real-time Quality Control (rQC®) system
Connection Interfaces	WLAN, USB, Ethernet interface, RS232 Serial Interface, compatible with LIS
Print Mode	Built-in thermal printer, office printer (optional), MDMP, LIS
Data Capacity	50,000 test results and quality control data
Operation Environment	Temperature 10-30°C, Humidity 40-85%

Celercare™ M5 Reagent Discs

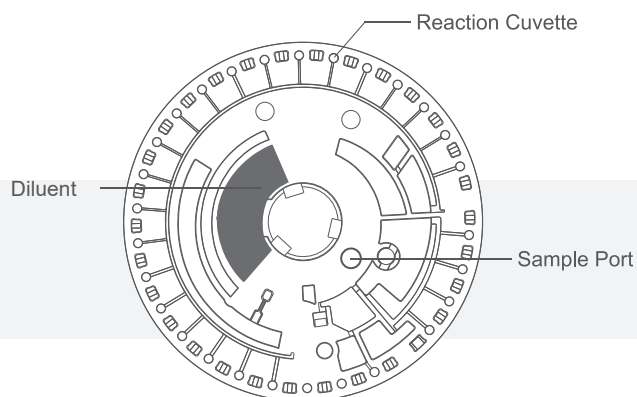

	General Chemistry I Lyophilized Kit	Clinical Emergency Lyophilized Kit	Renal Function Panel Lyophilized Kit	Liver Function Panel Lyophilized Kit	Electrolyte Panel Lyophilized Kit	Glucose and Lipid Panel Lyophilized Kit	General Chemistry II Lyophilized Kit	Liver and Renal Function Lyophilized Kit	Ammonia Panel Lyophilized Kit	General Chemistry IV Lyophilized Kit
TP	TP			TP				TP		TP
ALB	ALB		ALB	ALB				ALB		ALB
GLO*	GLO*			GLO*				GLO*		GLO*
A/G*	A/G*			A/G*				A/G*		A/G*
ALT	ALT			ALT				ALT		ALT
AST	AST	AST		AST				AST		AST
GGT				GGT				GGT		GGT
ALP				ALP						ALP
TBIL	TBIL			TBIL				TBIL		TBIL
DBIL	DBIL			DBIL						DBIL
IBIL*	IBIL*			IBIL*						IBIL*
UA	UA	UA	UA							UA
CRE	CRE	CRE	CRE				CRE	CRE		CRE
UREA	UREA		UREA				UREA	UREA		UREA
CK		CK								
CK-MB		CK-MB								
LDH		LDH								
α-HBDH		α-HBDH								
AMY		AMY					AMY			
GLU	GLU	GLU				GLU	GLU	GLU		GLU
GSP						GSP				
TG	TG					TG				TG
CHOL	CHOL					CHOL				CHOL
HDL-C	HDL-C					HDL-C				HDL-C
LDL-C*	LDL-C*					LDL-C*				LDL-C*
K ⁺		K ⁺			K ⁺		K ⁺			
Na ⁺		Na ⁺			Na ⁺		Na ⁺			
Cl ⁻		Cl ⁻			Cl ⁻		Cl ⁻			
Ca ²⁺			Ca ²⁺		Ca ²⁺					
P			P		P					
Mg ²⁺					Mg ²⁺					
CO ₂		CO ₂	CO ₂		CO ₂		CO ₂			
NH ₃									NH ₃	

* Calculated test value